

## Vital patient data at a glance



**Sentec Connectivity**  
enables the transmission of  
oxygenation and ventilation  
monitoring – **tcPCO<sub>2</sub>**, **tcPO<sub>2</sub>** –  
to Dräger Infinity®  
patient monitors.




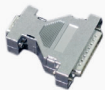

## Sentec Connectivity Dräger Infinity® patient monitors

- Sentec Digital Monitor readings are duplicated to Infinity patient monitors
- The Sentec Digital Monitor is automatically identified when connected to an Infinity patient monitor
- Cost-effective data integration

Care with  
Confidence

# Sentec Connectivity

Dräger Infinity® patient monitors

| Required Components   | System Requirements  | Description  |
|---|--|--|
| <b>Sentec Digital Monitor (SDM)</b><br>                            | Software version V08.00 to V08.03 (using the „TCB“ protocol) | Provides noninvasive and continuous real-time monitoring of ventilation and oxygenation  |
| <b>Infinity patient monitors</b><br>                               |  | Displays and stores tcPCO <sub>2</sub> and tcPO <sub>2</sub> on<br>Infinity Delta<br>Infinity Delta XL<br>Infinity Kappa<br>Infinity Omega serie |
| <b>Medical Information Bus Connector (MIB Connector)</b><br>       |  | The MIB Connector is located in the Dräger Docking Station   |
| <b>Medical Information Bus II Converter (MIB II Converter)</b><br> |  | Connects the patch cable with the Infinity Adapter   |
| <b>Infinity Adapter</b><br>                                       |  | Connects the SDM with the MIB II Converter   |
| <b>Patch cable</b><br>   |  | Connects the MIB II Converter and the MIB Connector to the Dräger Docking Station  |
| <b>Standard RS232 cable</b><br>                                  |  | Connection between SDM and Infinity Adapter  |

## Sentec Connectivity Support Team

Let's connect.

If you need support, please contact your local Sentec representative, send an e-mail to [connectivity@sentec.com](mailto:connectivity@sentec.com) or consult your Dräger representative.

Please visit [www.sentec.com](http://www.sentec.com) for more information.

Your local distributor:

**sentec.**

Ringstrasse 39  
4106 Therwil  
Switzerland  
[www.sentec.com](http://www.sentec.com)  
 Made in Switzerland

Sentec Connectivity Support Team:  
Tel +41 61 726 97 60  
[connectivity@sentec.com](mailto:connectivity@sentec.com)

**Dräger**

Drägerwerk AG & Co. KGaA  
Moislinger Allee 53-55  
23558 Lübeck  
Deutschland  
Tel +49 451 882 0  
[www.draeger.com](http://www.draeger.com)

**Note:** This is a global overview and may contain information on products not released in all markets/countries. Information is subject to change without notice.